



Ermis[®]

MedTech



Medical + Optical

Caring about bringing quality to you

info@medopt.com.au

Toll Free: 1800 225 065

www.medopt.com.au



Ermis[®]

MedTech

Manufactured in Germany since 1998 from Aluminium alloy with an anodised surface to protect against corrosion.

A great range of container heights available to enable use of a variety of tray heights, light weight and with a 2 year warranty.

Machine wash for reprocessing using a suitable neutral cleaning agent.

Plastic Security Seal with locking tag and steam indicator, must be attached to the closures outside, before each sterilisation.

The indicator changes its colour from pink to brown, during steam sterilisation at 134°.

ERMIS* sterilisation containers have been tested for a shelf life of 12 months.



E Standard, paper filter with steam indicator



Premium System 3-In-One Aluminium filter



TUT Premium System



Container Bottoms



Bottom Container 1/2
ER305.100E
297x274x88 mm



Bottom Container 1/1
ER105.150E
585x274x138 mm



Bottom Container 3/4
ER205.135E
462x274x123 mm

Container Lids

We offer a range of aluminium lids to suit your unique requirements and budget.

All lids are designed to fit the three sizes and various heights of container bases illustrated above and come in a variety of colours for easy specialty identification.

E Standard, paper filter with steam indicator

Premium System 3-In-One Aluminium filter

TUT Premium System

E Standard Model

Classic single use paper filter with steam indicator



Container Lids

Premium System 3- In-One Aluminium Filter

Barrier Filter, life long & maintenance free

Based on the isolation of airflow & micro-organisms




Container Lids

TUT Premium System

Sterile Container with aluminium solid lid and unique labyrinth barrier filter

Life long & maintenance free

Based on the principle of the Pasteur loop

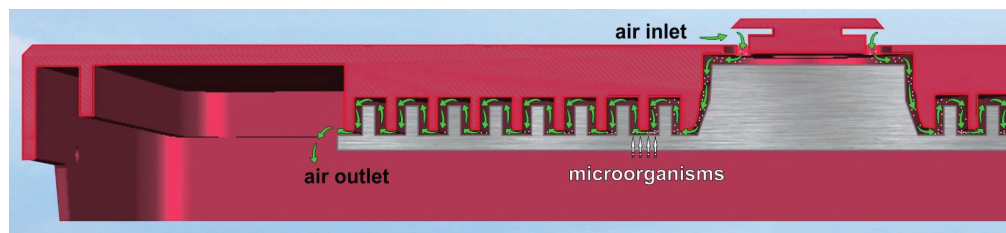

Accessories



Aluminium Tag
with option to
laser mark for
easy tray
identification

ER821.030B
Red tag




Plastic Single Use Security Seal with
locking tag and steam
indicator changing colour from pink
to brown during steam sterilisation at
34°.

ER818.010 Blue pack/100



The Silicone gaskets in the lids can be
changed easily and needs to be checked
regularly and changed after 500 cycles

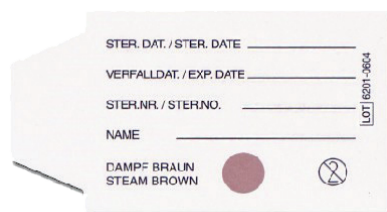

Single use paper
labels with steam indicator
ER810.000E Pack 100



Single use paper
filter with steam indicator
ER800.000E Pack 100



Product Index

Container Bottoms					
1/2 Bottom Container					
	ER305.100E	ER305.135E	ER305.150E	ER305.200E	ER305.260E
SIZE	297x274x88mm	297x274x123mm	297x274x138mm	297x274x188mm	297x274x248mm
3/4 Bottom Container					
	ER205.100E	ER205.135E	ER205.150E	ER205.200E	
SIZE	462x274x88mm	462x274x123mm	462x274x138mm	462x274x188mm	
1/1 Bottom Container					
	ER105.100E	ER105.135E	ER105.150E	ER105.200E	ER105.260E
SIZE	585x274x88mm	585x274x123mm	585x274x138mm	585x274x188mm	585x274x248mm

Container Lids			
E-Standard			
	1/2 Size	3/4 Size	1/1 Size
	ER300.010E	ER200.010E	ER100.010E
	ER300.020E	ER200.020E	ER100.020E
	ER300.030E	ER200.030E	ER100.030E
	ER300.040E	ER200.040E	ER100.040E
	ER300.050E	ER200.050E	ER100.050E
	ER300.060E	ER200.060E	ER100.060E
TUT Premium System			
	1/2 Size	3/4 Size	1/1 Size
	ER300.010TB	ER200.010TB	ER100.010TB
	ER300.020TB	ER200.020TB	ER100.020TB
	ER300.030TB	ER200.030TB	ER100.030TB
	ER300.040TB	ER200.040TB	ER100.040TB
	ER300.050TB	ER200.050TB	ER100.050TB
	ER300.060TB	ER200.060TB	ER100.060TB

Accessories			
Aluminium Coloured Identification Tags			
	No Marking	With Markings	
	ER821.010	ER821.010B	
	ER821.020	ER821.020B	
	ER821.030	ER821.030B	
	ER821.040	ER821.040B	
	ER821.050	ER821.050B	
	ER821.060	ER821.060B	
	ER821.070	ER821.070B	
	ER821.080	ER821.080B	
	ER821.100	ER821.100B	





Ermis[®]

MedTech



Mini & Maxi Containers

Medical + Optical

Caring about bringing quality to you

info@medopt.com.au

Toll Free: 1800 225 065

www.medopt.com.au

Mini & Maxi Containers



ER405.100E 310x131x85 mm



ER400.010E



ER400.010TB



ER505.060E 312x183x36 mm




ER500.010E



ER500.010TB



ER507.073E 227x161x73 mm

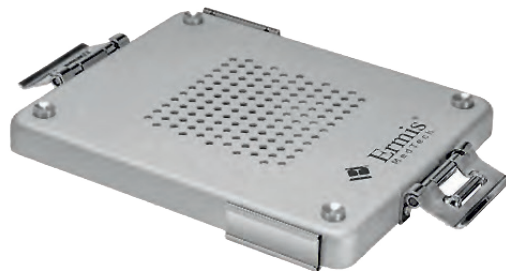

ER509.010E



ER509.010TB



ER505.50E 211x140x41 mm



ER504.010E



ER504.010TB




ER805.000



ER802.000E



ER818.010
Blue pack 100

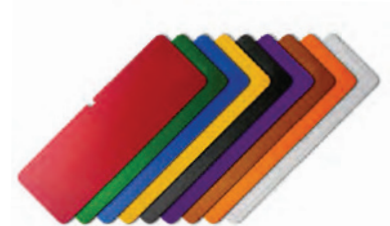

Aluminium Tag
ER820.010



Product Index

Container Base		
ER405.040E	ER405.070E	ER405.100E
310x131x28mm	310x131x56mm	310x131x85mm
Mini lid perforated paper filter		
	ER400.010E	
	ER400.020E	
	ER400.030E	
	ER400.040E	
	ER400.050E	
	ER400.060E	
TUT Premium lid, multi use labyrinth barrier filter		
	ER400.010TB	
	ER400.020TB	
	ER400.030TB	
	ER400.040TB	
	ER400.050TB	
	ER400.060TB	

Container Base		
ER505.050E		
211x140x41mm		
Mini lid perforated paper filter		
	ER504.010E	
	ER504.020E	
	ER504.030E	
	ER504.040E	
	ER504.050E	
	ER504.060E	
TUT Premium lid, multi use labyrinth barrier filter		
	ER504.010TB	
	ER504.020TB	
	ER504.030TB	
	ER504.040TB	
	ER504.050TB	
	ER504.060TB	

Container Base				
ER505.040E	ER505.060E	ER505.080E	ER505.100E	ER505.130E
312x183x36mm	312x183x54mm	312x183x72mm	312x183x85mm	312x183x122mm
Mini lid perforated paper filter				
	ER500.010E			
	ER500.020E			
	ER500.030E			
	ER500.040E			
	ER500.050E			
	ER500.060E			
TUT Premium lid, multi use labyrinth barrier filter				
	ER500.010TB			
	ER500.020TB			
	ER500.030TB			
	ER500.040TB			
	ER500.050TB			
	ER500.060TB			

Container Base			
ER507.060E	ER507.073E	ER507.100E	ER507.130E
230x161x55mm	230x161x73mm	230x161x90mm	230x161x120mm
	4 lg tag holders		
Mini lid perforated paper filter			
	ER509.010E		
	ER509.020E		
	ER509.030E		
	ER509.040E		
	ER509.050E		
	ER509.060E		
TUT Premium lid, multi use labyrinth barrier filter			
	ER509.010TB		
	ER509.020TB		
	ER509.030TB		
	ER509.040TB		
	ER509.050TB		
	ER509.060TB		

Accessories		
Aluminium Coloured Identification Tags		
	No Marking	With Markings
	ER820.010	ER820.010B
	ER820.020	ER820.020B
	ER820.030	ER820.030B
	ER820.040	ER820.040B
	ER820.050	ER820.050B
	ER820.060	ER820.060B
	ER820.070	ER820.070B
	ER820.080	ER820.080B
	ER820.100	ER820.100B
Blue Security Seal	ER 818.010	Pack 100
Paper filter square	ER805.000	Pack 100
Paper filter oblong	ER802.000E	Pack 100



Hard to Find Instruments?

Let us do the looking for you



Medical + Optical Head Office
81 Whiting St
Artarmon NSW 2064
info@medopt.com.au
TOLL FREE: 1800 225 065
www.medopt.com.au
ABN: 80 003 178 596



Ermis[®]
MedTech



KLS martin
GROUP

Surgistar[™]
The Cutting Edge of Tomorrow

